

CuroCell Cirrus®
Alternating pressure system

User manual item number: 95-001054



Contents

Important safety precautions	3
Introduction	4
Product description	8
Assembly/Installation	11
Operation	12
Cleaning and reconditioning	14
Storage	15
Maintenance	16
Troubleshooting	17
Technical description	18
Spare parts	20
Options and accessories	21
Warranty	22
Other information	23

IMPORTANT SAFETY PRECAUTIONS

Read all instructions prior to use

WARNING! To minimise the risk of fire, electric shocks or bodily injury:

- 1 The product must only be used for its intended purpose as described in this manual. The product may not be combined, assembled or repaired with other parts (e.g. pump and mattress), accessories or spare parts other than those described in this manual or other documentation from Care of Sweden.
- 2 Do not use the pump in wet rooms or where there is a risk of the pump coming into contact with water or another liquid. Never handle a pump that has come into contact with water/liquid. Pull the plug out of the electrical socket immediately and send the product to an authorised service technician for servicing.
- 3 DO NOT USE in the vicinity of or in contact with fire sources/hot surfaces, such as fire, burning cigarettes or hot lamps. Although the product has undergone fire testing, it could become damaged if it comes into contact with fire sources.
- 4 Position and use the product in a manner that prevents it from becoming trapped or damaged. NOTE: Be particularly aware of trapping damage when using bed rails.
- 5 Never use the product if the pump power cord or plug is defective, the pump casing is damaged or it is not functioning properly. Contact a service technician for examination and repair.
- 6 Regularly check product functionality through a hand check (see section 4.4). Also check that the comfort setting is correct as indicated in the table on the pump.
- 7 Continual monitoring is required when the product is used for individuals needing special supervision, such as children.
- 8 The mattress are protected with a hygiene cover; Please try to avoid the simultaneous use of multiple hygiene cover with this mattress. Because of the risk vapour permeable can be destroyed.
- 9 The hygiene cover does not allow liquid or air to penetrate, but is vapour permeable. Make sure that the patient is positioned correctly to avoid the risk of suffocation.
- 10 Prevent damage to the hygiene cover from sharp objects.
- 11 If the cover is equipped with side handles, these are for lifting the mattress. All other use takes place under your own liability and is not covered by the product guarantee.

1. Introduction

CuroCell® air mattresses operate according to the alternating pressure principle, a technique that is well proven and been used for a long time. Alternating pressure is when the pressure regularly alternates between the mattress cells to replicate body movement. Variable pressure cycles are set up where there is little or no pressure. This is optimal for skin improvement or to prevent skin damage.



Always read the user manual prior to use.

1.1 General information

This CuroCell® mattress is a medical device in accordance to MDD 93/42/EEC. According to statutory regulations made by the authorities regarding medical devices, the manufacturer is required to report all accidents or incidents involving there products. We would be very grateful for all information involving accidents or incidents relating to our products, by reporting them immediately to us, here at Care of Sweden. The system has been tested and approved according to the following standards: IEC 600601-1, IEC 60601-1-2, SS-EN ISO 12182, SS-EN ISO 14971, SS-EN ISO 10993, EN 597-1, EN 597-2, EN 14126.

1.2 Intended use

CuroCell Cirrus® combines alternating pressure with foam and can be used in all types of health environments, to prevent pressure ulcers, and as an aid in the treatment of category 4⁽¹⁾(*). It is suitable for lightweight and sensitive people or to people with pain. The mattress has a built in heel function for added relief of heels. The mattress is placed directly on the bed base. The recommended user's weight is up to 230 kg.

For the Evac-function the maximum user weight is 130 kg.



Important! (*) For users between 180-230 kg the recommended use is changed to treatment of pressure ulcers up to category 2⁽¹⁾.

Cleaning of the mattress shall be made between replacements of each user.

We recommend that position changes are made regularly. The time interval must be evaluated by the responsible personnel, depending on the patient's status diagnosis and general condition.

NOTE!

- When used together with positioning pillows, read the pillow user manual carefully for correct positioning of the patient.
- If you are using this product with an evacuation sheet, it is the user's responsibility to ensure that it is safe to use for evacuation.
- Be careful when using bed rails or other protection on the bed so that the mattress does not become trapped or damaged.
- The mattress must not be lifted with a person lying on it or used for any other type of transport other than a possible evacuation using an evacuation sheet.
- This mattress is unsuitable for x-ray examinations because of the risk of blurred images or artefacts that may lead to diagnostics error.

(1) National Pressure Ulcer Advisory Panel, European Pressure Ulcer Advisory Panel and Pan Pacific Pressure Injury Alliance. Prevention and Treatment of Pressure Ulcers: Quick Reference Guide. Emily Haesler (Ed.). Cambridge Media: Osborne Park, Western Australia; 2014.

1.3 Other information

Action plan Pressure ulcers

To achieve the best possible result for this product it is recommended that a well-planned and well-structured working environment be used. In our guideline "Action plan" you will find examples which can be applied and put into practice when using our products.

Mattress guide

This practical guide gives you as a care giver, help in choosing the correct mattress from the Care of Swedens range. CuroCell Cirrus® is a mattress in function group C. This information can be found on our website www.careofsweden.com or you can order a copy from our Customer services by sending your details to info@careofsweden.se.

Evac-cover

A cover with protected zipper on all sides, divisible (top and bottom part).
Functional handles for evacuation, moving and transport. Colour grey/black.



1.4 Evacuation instructions


Evac is a mattress cover with a built-in function for evacuation of bed-bound users. In case of fire or other emergency, there is no time to think about what to do next. Therefore, it is very important that the safety equipment, as well as the evacuation procedures, are well planned and known in advance.

The Evac cover is equipped with three handles on the head and foot end. On the sides, there are handles with Velcro straps. These straps are to be used to fixate the user to the mattress before commencing the evacuation. The straps are folded in the protective pockets on the side of the mattress with the ends sticking out to ensure that they are easy to pull out.

WARNING!

Person to perform an evacuation needs to have the necessary training. Always check the equipment, cover and mattress before commencing, to make sure an evacuation can be carried out safely. It is the user's responsibility to ensure that it is safe to use for evacuation.

Evacuation instructions

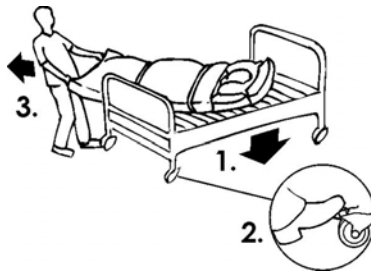

Down from the bed

- Inform the patient about what will happen.
- Lower and lock the bed in position.
- Disconnect the tubes from the pump.
- Pull out and fasten the Velcro straps.
- Use the outer handles at the foot end of the mattress and turn the mattress out from the bed.
- Slowly pull the patient down on to the floor.



On an even surface

- Pull on the outer handles at the foot end of the mattress.
- Working position should be a little bit leaning back.

Down the stairs

- Walk backwards until the main part of the mattress has slide pass the top stair.
- Turn the back to the patient and look in the moving direction.
- Hold one hand on the stair handle and pull with the other hand on the middle handle of the mattress.
- Walk down the stairs. It is important to maintain speed, without running.



Training/practise

To bear in mind when training.

Pulling down from the bed:

- The patient in the bed should not be at the very top of the mattress. There could be a risk of hitting the head in the bed when pulling down.

Pulling in the stairs:

- Start with pulling an empty mattress down the stairs, in order to practice the technique.
- When training with a person on the mattress, a string or rope can be fasten in the handles at the head end. Another person can then walk behind the mattress and help out if problems should occur.

- It is advisable to practice in the same stairs where an evacuation can be necessary.

General:

- The main purpose of the training is not to reach high speed, but to learn the correct technique.

Training covers

- Because of the heavy strain on the equipment when evacuating, always check the seams and handles at the end of training. We recommend that equipment used for training is marked with “Training equipment” and is only used for this purpose.

Explanations

- If the Velcro straps are not used, there is a risk of the patient sliding/rolling down from the mattress.
- Pulling at the foot end reduces the strain on the puller, which means faster and safer evacuation is possible.
- Pulling in the stairs is the most difficult part of the evacuation and should be practiced extra carefully.
- The arms of the patient shall be inside of the Velcro straps. This way reduces the risk of the user grabbing hold to something out of fear.

NOTE!

There is a difference in maximum user weight between use of the mattress and use of the Evac function. Please read section 1.2 for more information.


To ensure proper function of the Evac cover, it is important to regularly check the quality of the mattress. Check for damages that could cause problems when evacuating.



2. Product description


2.1 Pump and mattress

CuroCell Cirrus® has eight alternating air cells and these cells are encased in HR-foam (high resilient foam). On top of the HR-foam, closest to the user, is a thin layer of visco-elastic foam.

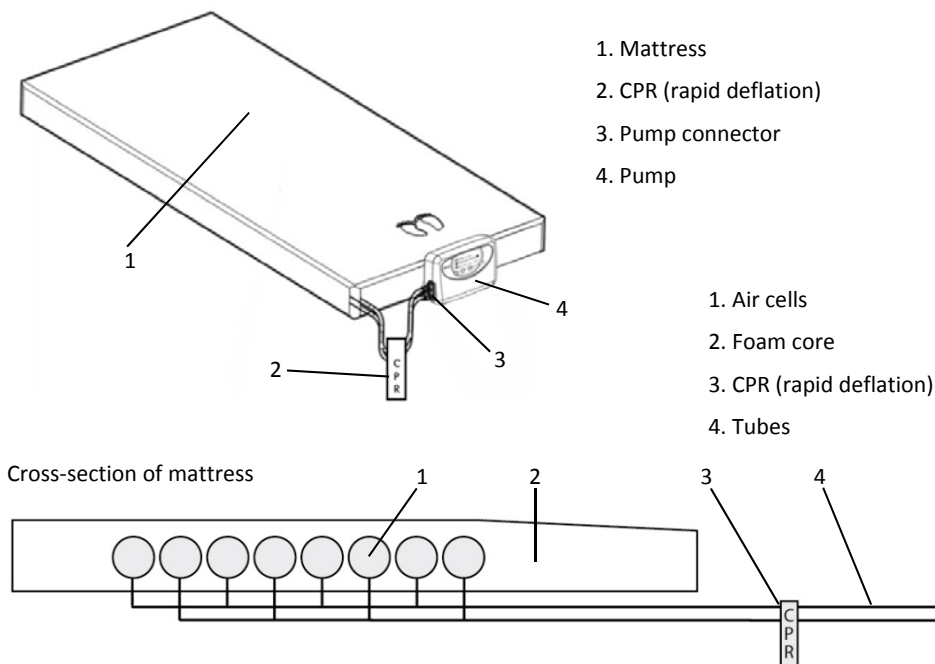


2.2 Pump

The CuroCell Cirrus® pump applies alternating pressure that is adapted to the weight of the user. The air cells inflate and deflate with different cycle times.

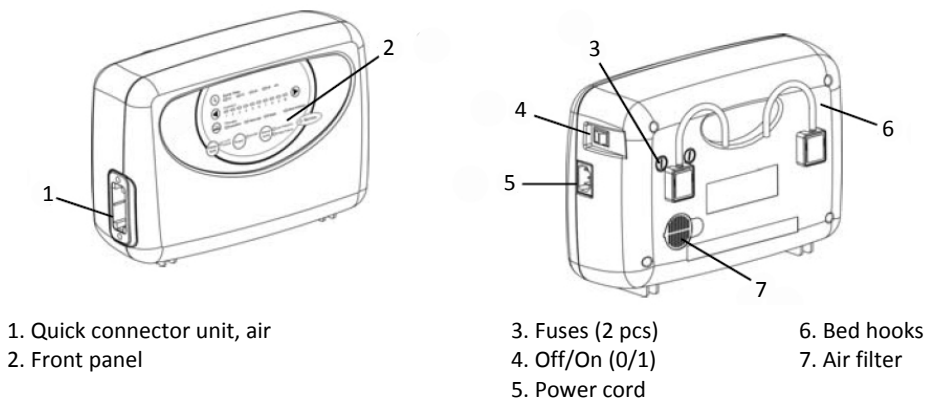

2.3 Pump functions – control panel

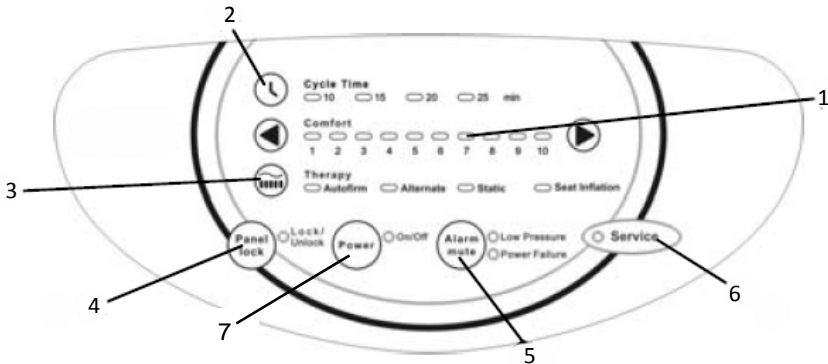

1 COMFORT SETTING

Suitable air pressure is set using the comfort dial. Guidelines for user weight can be found in the table on the pump.

2 CYCLE TIME

There are four selectable cycle times. Press the button to select cycle time.

3 PROGRAM (THERAPY)

Change the mode by pressing the button, the chosen mode is indicated by the green light. The modes' order is: Autofirm, Alternate, Static, Alternate - seat inflation, Static - seat inflation.

A. Autofirm Autofirm is used to fill the mattress with air completely for firm support. The pump always starts in this mode. When the correct pressure is reached, it changes automatically to the previous setting or to the pre-setting outlined below.

Presetting features for Autofirm process

1. Program (Therapy): Select the required program.
2. Comfort setting: Set the required air pressure (1-10).
3. Cycle time: Select from 10, 15, 20 or 25 minutes.

Note: The preset Cycle Time, Comfort Level and Therapy mode will be illuminated to indicate the settings during the Autofirm process. (Autofirm LED is also illuminated.) After mattress is inflated, Autofirm LED will be turned off automatically.

B. Alternate In this mode, the pressure alternates between the cells and there are four different cycle times to choose from.

C. Static When "Static" is chosen, the alternating process is terminated and the pressure between the cells evens out. The function can be used in association with patient care, repositioning/transferring in and out of bed. The function is turned off by pressing the button again, or it will automatically go back to alternating mode after 20 minutes.

D. Seat inflation This mode gives extra support when sitting up and can be set for both alternating and static mode.

4 PANEL LOCK

Press down the Panel lock button for about 2 seconds to lock or unlock the panel.

5 ALARM MUTE

Press the ALARM MUTE button to temporary turn off the alarm. Should the situation not be resolved within 5 minutes, the alarm will be activated again. Please refer to section 8 to troubleshoot any problems with the alarm.

A. Power failure alarm

During a power failure, both the sound and light alarm will go on. Turn it off by pressing the ALARM MUTE button.

B. Low pressure alarm

When the low pressure light is on, the pressure in the mattress is below normal. Both the light and sound alarm will activate. The alarm will sound until the ALARM MUTE button is pressed. If the problem remains, see chapter 8 for trouble shooting.

6 SERVICE


This light indicates a technical failure in the pump. First check that the air filter is clean and that the pump is not covered by anything (such as a pillow or a blanket). If the problem remains, please contact the manufacturer or your local dealer for further information.

7 POWER

Power on/off.

3. Assembly/Installation


When unpacking the mattress, check that none of the parts are damaged. In the event of damage, contact the manufacturer or your local dealer before using the product.

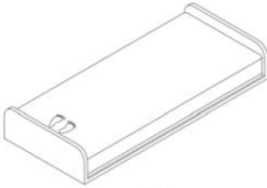

1. Place the mattress on the bed base.


Note!

- Make the bed with sheets for added comfort.
- Check that the CPR valve is closed.

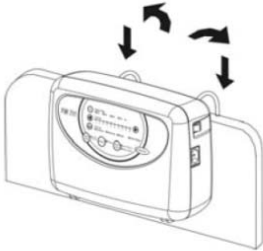

2. Hang the pump at the foot end of the bed or place it on a flat, sturdy surface.

Route the power cord so that there is no risk of it causing someone to trip, being run over by the bed wheels or getting pinched when the bed is raised or lowered.

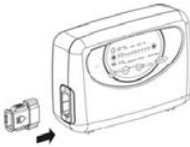

3. Connect the air tubes from the mattress to the side of the pump. There is an audible “click” sound when the quick connectors are properly attached.



4. Check that the switch on the pump is in the “0” position (Off). Connect the plug to an approved electrical socket (220 V).

4. Operation



Note! Always read the user manual prior to use.

4.1 Start-up

1. Set the power switch on the side of the pump to 1 (On) and start the pump by pressing the power button on the front panel.
2. Air starts to fill up the mattress. It takes about 20 minutes until the mattress is ready for use. **Note!** The pump always starts in Autofirm mode.
3. The light for 'Low pressure' is lit until the mattress has reached the set pressure. When the light goes out, the system changes mode (to previous settings) and the mattress are ready to be used.



4.2 Program selection

There are four different programs. Please refer to section 2.3 for more details. As a standard setting, we recommend an alternating mode and a 10 minute cycle time.

4.3 Comfort settings

The air pressure in the mattress can be adjusted according to the patient's desire for comfort. ALWAYS consult the prescriber or other authorised personnel for the appropriate setting.

The mattress is now ready for use. Use the comfort buttons (1-10) to set a suitable air pressure based on the table on the pump.

Note! The pump always starts in the Autofirm mode. When the mattress has reached the required pressure, the pump is set automatically to the previously used setting or to the settings made in the Autofirm program. Please refer to section 2.3 for more details.

kg	
0-40	1
50	2
60	3
70	4
80-90	5
100-110	6
120-140	7
150-160	8
170-180	9
190-230	10

NOTE!

- We recommend the Seat inflation mode for all situations where the head of the bed is raised.
- When only parts of the mattress are under load, e.g. if the patient is an amputee, the pressure should be increased (1-2 levels).
- Always check the set pressure with a "hand check", see 4.4 for more information.

4.4 Hand check

Hand check is used to ensure that the mattress is properly inflated. This should be done regularly, once per shift and after changing the position in the bed or changing the comfort settings.

Insert your hand between the air cells directly underneath the user's pelvic region, use a vertical hand. There should be a clearance of 2-3 fingers between mattress base and the patient. If not, you need to adjust the pressure setting.



4.5 CPR – Rapid deflation of the air cells

To empty the air cells, pull the CPR (Cardio pulmonary rescue) band placed between mattress and pump. In case of an emergency where CPR is necessary, follow the local instructions.

Note! Make sure that the CPR connector is properly connected after use.



4.6 Alarms

When the pressure drops to an abnormally low level, the Low Pressure diode lights up together with an audio alarm. Please refer to section 8 for troubleshooting. The alarm is deactivated by pressing the ALARM MUTE button on the pump's control panel. If the problem is not resolved within 5 minutes, the alarm will start again. The diode remains on until the pressure level is remedied and reset.

The pump is also fitted with a Power Failure alarm. An alarm sounds and flashes when the pump is disconnected from the electric power. Turn off the alarm by pressing the 'ALARM MUTE' button.

WARNING!

If the air pressure remains low, check the mattress and tubes for damage. Repair if necessary or contact the manufacturer or your local dealer for help.

4.7 Transport function

If the patient must be moved in bed, the tubes can be connected together (detached from the pump) and the mattress put in "transport mode". We recommend that this function is only used for short periods of time.



5. Cleaning and reconditioning

It is important to follow the instructions below before using the product for another user. Consult the Clinician responsible or the supplier for help and instructions when unsure.

5.1 Cleaning

WARNING!

Do not open the pump housing – risk of electric shock. Service and maintenance must be performed by authorised technicians. Do not get the pump wet or submerge it in any liquid.

PUMP



Wipe of the pump with a damp cloth and a mild cleaning solution for this purpose such as: alcohol with or without tensides, Isopropanol 70 %, oxidating solutions or chlorine max 1 %. If using a different agent, choose one that does not affect the pump housing.

MATTRESS COVER



Covers can be wiped off with a mild cleaning solution for this purpose such as: alcohol with or without tensides, Isopropanol 70 %, oxidating solutions or chlorine max 1 %.

Excrements and blood spots must be removed immediately with cold water. Carefully follow local instructions and usage instructions for the cleaner.



Multi-layer covers should be divided before washing. Wash with like colours.



Chlorine and phenol-based cleaners could have a negative impact on the PU surface and should be avoided. If chlorine is used, we recommend a mixture of max. 1%.

INNER COVER SOFT



Clean the affected area with a mild cleaning solution for this purpose such as: alcohol with or without tensides, Isopropanol 70 %, oxidating solutions or chlorine max 1 %.

5.2 Reconditioning

PUMP

See section 5.1 for instructions.

MATTRESS

1. Remove the covers. See instructions in chapter 5.1 for cleaning covers.
2. Disconnect the tubes from the pump.
3. Carefully pull the air cells out from the foam core.
4. Place the air cells and the tubes on a clean area which is suitable for washing.
5. Wipe off the cells, all tubes and the CPR valve with a mild cleaning solution for this purpose such as: alcohol with or without tensides, Isopropanol 70 %, oxidating solutions or chlorine max 1 %.
6. Place the air cells back into the foam core. Check that all the tubes and connectors are secure and tight, before closing the cover.
7. Wipe off work surfaces with a mild cleaning solution for this purpose such as: alcohol with or without tensides, Isopropanol 70 %, oxidating solutions or chlorine max 1 %. If cells for any reason have loosened from the tubes, these must be replaced according to drawing in chapter 2.1.

FOAM CORE



Clean the affect area with a mild cleaning solution for this purpose such as: alcohol with or without tensides, Isopropanol 70 %, oxidating solutions or chlorine max 1 %.

Press out the water carefully. Let the foam core dry in a warm, ventilated area, without direct exposure of sunlight. The foam core needs to be completely dry before it is used again.

NOTE!

Check the hygiene cover each time the product is cleaned. If damaged, it must be replaced or repaired. Also check the pump, tube connectors and power cord during cleaning. Damaged parts must be replaced or repaired.

6. Storage


Store the mattress on a flat surface or carefully rolled/folded with the zipper closed. Always store at room temperature in a dry area. Make sure the mattress is stored where it will not be exposed to thrusts, impact or other things that can damage it. Avoid direct exposure of sunlight. If the mattress is to be stored more than 3 months, we recommend it is stored on a flat surface (not folded or rolled up).

7. Maintenance

7.1 General


1. Check the power cord and plug for damage and scratches.
2. Check that the mattress cover is intact and that the cover and air cells are correctly assembled.
3. Check that the air flow alternates between each outlet on the side of the pump.
4. Check for damaged or pinched tubes and connectors. Contact the manufacturer or your local dealer for spare parts.

7.2 Changing fuses

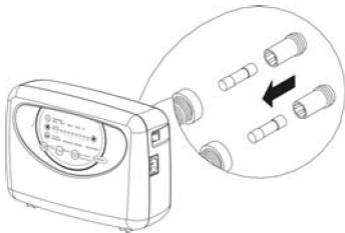

1. Pull the plug out of the electrical socket before replacing fuses.
2. Remove the cover from the fuses with a screwdriver. Pull out the fuse holder.
3. Fit a new fuse of the same type (T1AL, 250V) and screw the same type and screw the protective cover back in place. There are two fuses.

7.3 Changing the air filter

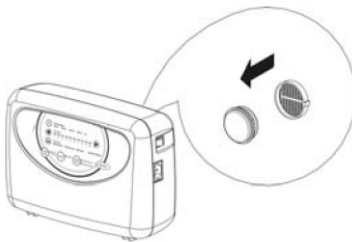

1. Unscrew the round safety plate on the back of the pump.
2. Lean the pump backwards in order for the filter to fall out from the holder.
3. Place the new filter in the filter holder and replace the safety plate.

If the pump is used in a dirty environment, the filter should be checked regularly.

8. Troubleshooting

Problem	Solution
Pump does not start	<p>Check that the power cord is connected to an approved electrical socket.</p> <p>Check the fuses.</p> <p>Make sure the power switch on the pump side is set to "1".</p>
Alarm is on and mattress is losing air	<p>Check that the quick connectors on the air tubes are properly connected to the side of the pump.</p> <p>Check that the CPR valve is closed.</p> <p>Check that all tube connectors are properly seated and connected.</p> <p>Open the mattress and check for leaks from air cells or tubes.</p>
The patient is bottoming out	<p>The pressure might be incorrectly set. Adjust according to the weight table and wait a few minutes. Check again using a Hand Check (see section 4.4).</p> <p>Check that the CPR valve is closed.</p>
Some cells have less air	<p>This is normal since the cells are alternately inflated and deflated when the system is set to alternating mode.</p>
Pump is noisy, vibrations are felt	<p>Check that the pump hangs correctly on the bed. Resonance can occur in parts of the bed. Remove the pump and listen whether this makes a difference. Can be resolved by putting the pump on a flat, steady surface or by placing a towel between the pump and bed.</p>
If your questions are not answered by the above information, please contact Care of Sweden or your local dealer.	

9. Technical description

9.1 Technical specification

PUMP		SPECIFICATION
Model		CuroCell® Vari
Voltage		220 V-240 V, 50-60 Hz
Current consumption		0.07 A
Fuse type		T1AL, 250 V
Cycle time		10, 15, 20 or 25 min
Dimensions (BxLxH)		29.1 x 11.7 x 20 cm
Weight		2.2 kg
Surroundings	Temperature	Operations: 10°C - 40°C
		Storage: -15°C - 50°C
		Shipping: -15°C - 70°C
	Humidity	Operations: 10% - 90% non-condensing
		Storage: 10% - 90% non-condensing
	Athmospheric pressure	700hPa-1060hPa
Classification		Class II, Type BF, IP41
MATTRESS		SPECIFICATION
Model		CuroCell Cirrus®
Dimensions (BxLxH)		80/85/90/105/120x200x20 cm
Weight		14 kg (width 85 cm)

NOTE: The manufacturer retains the right to modify the product specification at any time.

WARNING!

Follow the operations and maintenance instructions carefully in order to prevent bodily injuries, as well as damages to the product or other equipment.

9.2 Labelling

The mattress are labelled with product name, production date, article number, size, care and wash instructions, EAN-bar code, as well as space for own labelling. The product also displays which side is up/down and head/foot.

9.3 Symbol Definitions



Article No



Manufacturer

XXXX-XX-XX

Year-Month-Day



Foot placement



Read the user manual



Group of function



CE marked acc. to
MDD 93/42/EEC



Recycling



Not in domestic refuse,
follow recycling instructions



Recommended user weight



Patient information
– category



Directly on bed base



Recommended user
weight when using Evac
function



Counteract shear force



Heel function



Do not wash



Machine wash at 70 °C



Machine wash at 95 °C



Chlorine



Tumble dry



Do not tumble dry



Drip dry



Do not iron



Do not dry clean



Wipe off



Class II equipment
(Double Insulated)




Type BF



IP-class
(Ingress Protection Rating)



Temperature range

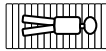

The mattress should only be
used lying in the lengthwise
direction on the mattress

10. Spare parts

Art. no.	Description	Size
CC-925531	Top part Steel	80x200 cm
CC-925532	Top part Steel	80x210 cm
CC-925533	Top part Steel	85x200 cm
CC-925544	Top part Steel	85x210 cm
CC-925535	Top part Steel	90x200 cm
CC-925541	Top part Steel	90x210 cm
CC-925542	Top part Steel	105x200 cm
CC-925549	Top part Steel	105x210 cm
CC-925539	Top part Steel	120x200 cm
CC-925529	Top part Steel	120x210 cm

CC-901001	Air cell, 80 cm wide mattress
CC-901002	Air cell, 85 cm wide mattress
CC-901003	Air cell, 90 cm wide mattress
CC-901004	Air cell, 105 cm wide mattress
CC-901005	Air cell, 120 cm wide mattress

CC-8217	Fuse (set of 2)	-
CC-8202	Air filter	-

If other spare parts are needed, please contact Care of Sweden or your local dealer for further assistance.

11. Options and accessories

11.1 Options

Sturdy transport bag

Transport and storage bag.

Top part CuroCell LF (low friction)

Top part with increased sliding function, can be used for example heavy, immobile patients.

11.2 Accessories

Cover Steel

Size	Bottom part Steel	Complete cover Steel	Inner cover Soft
80x200 cm	CC-925111	CC-925631	CC-928011
80x210 cm	CC-925112	CC-925632	CC-928012
85x200 cm	CC-925114	CC-925633	CC-928013
85x210 cm	CC-925115	CC-925644	CC-928010
90x200 cm	CC-925117	CC-925635	CC-928015
90x210 cm	CC-925118	CC-925641	CC-928020
105x200 cm	CC-925120	CC-925642	CC-928018
105x210 cm	CC-925131	CC-925649	CC-928029
120x200 cm	CC-925122	CC-925639	CC-928010
120x210 cm	CC-925133	CC-925651	CC-928031

Steel/Evac cover

Size	Top part, Steel	Bottom part Evac	Complete Steel/Evac
80x200 cm	CC-926041	CC-926210	CC-925841
80x210 cm	CC-926042	CC-926211	CC-925842
85x200 cm	CC-926043	CC-926212	CC-925843
85x210 cm	CC-926044	CC-926213	CC-925844
90x200 cm	CC-926045	CC-926214	CC-925845
90x210 cm	CC-926046	CC-926215	CC-925846
105x200 cm	CC-926047	CC-926216	CC-925847
105x210 cm	CC-926048	CC-926217	CC-925848
120x200 cm	CC-926049	CC-926218	CC-925849
120x210 cm	CC-926050	CC-926219	CC-925850

Top part CuroCell LF (low friction)

Art.no.	Size	Art.no.	Size
CC-916210	80x200 cm	CC-916215	90x210 cm
CC-916211	80x210 cm	CC-916216	105x200 cm
CC-916212	85x200 cm	CC-916217	105x210 cm
CC-916213	85x210 cm	CC-916218	120x200 cm
CC-916214	90x200 cm	CC-916219	120x210 cm

Art.no.	Description	Size
58-900905	Transport bag	105-120 cm
58-900906	Transport bag, extra high	80-90 cm
58-900911	Transport bag, extra high	105-120 cm

11.3 Other accessories

To ensure the best level of comfort, we recommend use of accessories from Care of Sweden since our products are designed to work well together. Please contact Care of Sweden or your local dealer for accessories or more information on our products.

12. Warranty

12.1 Scope

This CuroCell® system is covered by a 2-year warranty for manufacturing defects. The warranty does not apply to normal wear and tear or product impairment or damage stemming from negligence or improper handling/care.

12.2 Service and repairs

This CuroCell® mattress is a medical device, in accordance with MDD 93/42/EEC, that is used as an aid in preventing and treating pressure ulcers. For this reason, we recommend that the pump be regularly serviced and inspected to maintain functionality and performance. Like other technical devices, a pump that is properly cared for works better and lasts longer.

Service and maintenance must always be performed by Care of Sweden or one of its authorised technicians. Only use spare parts that have been approved by Care of Sweden. Use of non-approved spare parts invalidates the warranty.

A warranty claim must be submitted to Care of Sweden before the product is returned.

Please contact the manufacturer or your local dealer if you have questions regarding warranty and service.

12.3 Life expectancy

This product's expected life is estimated to be 5 years.

13. Other information

13.1 Disassembly and recycling

Except for certain parts of the pumps, energy recovery is possible for almost all material in CuroCell® products through incineration in waste incineration facilities.

Pump:

A used CuroCell® pump must not be disassembled. It must instead be taken to a recycling centre and sorted as electronic waste.

Mattress:

A used CuroCell® mattress must be taken to a recycling centre and sorted as combustible waste.

For more information contact Care of Sweden or your local dealer.

13.2 Returns and complaints

Complaints are to be reported to the supplier. Return of a product shall follow the supplier's instructions and can only be done after the complaint has been reported. Damages due to the transportation should be reported to the transportation agent.



Manufactured by

care of sweden
Supporting Life

Contact:

Phone: + 46 771 106600

Fax: + 46 325 12840

Email: info@careofsweden.se

Internet: www.careofsweden.com

Address:

Care of Sweden AB
P.O. Box 146
SE-514 23 Tranemo
SWEDEN

Visit address:

Fabriksgatan 5A
SE-514 33 Tranemo
SWEDEN

Cargo address:

Byns väg 4A
SE-514 33 Tranemo
SWEDEN

Supplied by